



Shown above: Wilson Frame Pad Set on Action[®] O.R. Performance Table Pads with Action[®] Prone Headrest, Composite Contoured Armboard, Small Overlay, and Chest Roll



Wilson Frame Pad Set

Designed to cover the Wilson Frame and protect patients from shear and pressure injury. These 1/2" thick pads, and robust 1" thick pads easily slip over frame and should be secured with Hook and Loop attachment.

For orthopedic and neuro procedures when Wilson Frame is used. Two-piece set.

40701	1/2" Thick	29 L x 10 W x 1/2" H	74 x 25.4 x 1.3 cm
40702	1" Thick	29 L x 10 W x 1" H	74 x 25.4 x 2.5 cm



Kambin Spinal Frame Pad Set

Designed to fit the original Kambin Spinal Frame.

This two-piece pad set fits over frame's patient supports. The 1/2" thick pads are secured to the frame by included hook & loop attachments positioned along the bottom of each pad as well as on tabs (2 each per pad) that fold over and attach to the outer edges of the frame.

40712	33 L x 13 3/4 W x 1/2" H	13 L x 5.41 x 1.3 cm
--------------	--------------------------	----------------------



Peg Board Pad Set

The Action[®] Peg Board Pad set complements most standard peg board positioning systems. This multi-piece set is designed to protect patient skin and soft tissue from damaging pressure and shear forces while in the high-pressure lateral position.

The medium-length overlay (Item #40100) lies atop the peg board (beneath the patient) while the adjustable peg post pads of varying lengths easily wrap around the pegs and are secured by hook and loop attachment. Overlay and peg pads are available as a set, or as individual units.

Set (#40719) includes:

- 1 Medium Overlay Pad, #40100
- 4 Short Peg Post Pads, #40719A
- 4 Tall Peg Post Pads, #40719B

Peg board and peg posts not included

40719	Peg Board and Post Pad Set	Set includes (1) 40100, (4) 40719A, (4) 40719B	
40719A	Short Post Pad	5 L x 5 W x 5 1/2" H	12.7 x 12.7 x 13.97 cm
40719B	Tall Post Pad	5 L x 5 W x 8 1/2" H	12.7 x 12.7 x 20.32 cm
40100	Medium Overlay	46 L x 20 W x 1/2" H	117 x 50 x 1.3 cm



Vac-Pac[®] Covers

This pad is designed to line the interior of the Olympic Vac-Pac[®] (bean bag) positioner to help mitigate patient pressure and shear while in the lateral position. It is applicable for orthopedic, thoracic or urology procedures that dictate the patient be stable in the lateral position.

Vac-Pac[®] is a registered trademark of Natus Medical Incorporated

40612	Vac-Pac [®] Cover Medium Size 30	37 L x 28 W x 1/2" H	93.98 x 71.12 x 1.27 cm
40612-31	Vac-Pac [®] Cover Large Size 31	35 1/2 L x 38 3/4 W x 1/2" H	90.17 x 98.43 x 1.27 cm
40613	Vac-Pac [®] Cover T-Shape Size 35	15 L x 34 W x 1/2" H (bottom) 9 L x 23 W x 1/2" H (top)	38.1 x 86.36 x 1.27 cm (bottom) 22.86 x 58.42 x 1.27 cm (top)

Fracture Table Hip Pad

Designed for use on the fracture table plate, the primary function of this pad is to reduce pressure and shear at the sacral area.



40703 Hip Pad



Shown with Post Pad

Fracture Table Post Pad

This 1/2" thick wrap-around pad reduces pressure and shear on the patient's scrotum and upper leg areas. It also eliminates the tedious chore of wrapping and rewrapping the fracture table post with various makeshift items such as gauze and towels.



40704 Post Pad

Fracture Table Hip Pad with Post Hole

This oval shaped hip pad features an internal 1.75" hole to accommodate the perineal post. Pad is designed to protect the patient's upper leg and scrotum areas from skin injury due to pressure or shear forces during use.



40707 Hip Pad w/ Post Hole



Shown with Post Pad

40703	Hip Pad	8 L x 5 1/2 W x 1/2 H"	20 x 14 x 1.3 cm
40704	Post Pad	7 L x 5 W x 1/2 H"	18 x 12.7 x 1.3 cm
40707	Hip Pad w/ Post Hole	3 1/8 x 3 7/16 x 1 1/8 OD"; 2 1/8 x 2 5/16 x 1 5/16 ID" 1 3/4" diameter hole to accommodate perineal posts	4.5 cm diameter hole to accommodate perineal posts



40705 Knee Pad

40706 Iliac Strap

Andrews Frame™ Knee Pad


Knee Pad designed for use with the Andrews Frame® on the knee platform to alleviate pressure and shear during use.

Andrews Frame™ Iliac Strap

Two-piece adjustable strap system featuring durable nylon fabric belt with removable polymer pad insert. Removable polymer pad reduces pressure and shear forces on iliac crest and upper thigh areas.

Polymer pad covers up to 24" across the hips. Strap is secured with adjustable hook and loop closure.

40705	Knee Pad	23 L x 16 1/4 W x 1 H"	56 x 40.6 x 2.5 cm
40706	Iliac Strap	24 L x 6 W x 3/8 H"	61 x 15 x .95 cm

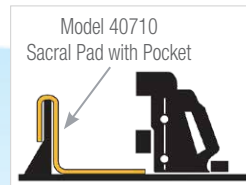

McGuire Hip/Sacral Pad

This 1/2" polymer pad protects the ilium, sacrum and area over the greater trochanter from shear and pressure damage while patient is in lateral position on the frame.

A pocket is created by folding the curved end of the pad over itself and securing attached hook and loop. To place on frame, slide pocket over the back plate, allowing rest of pad to cover the interior of the back plate and base of the frame.

40710	26 L x 10 W x 1/2 H"	66 x 22.5 x 1.3cm
--------------	----------------------	-------------------



Montreal Frame Pad

A positioning pad for the main platform of the Montreal Lateral Positioning Device. This pad is used in a lateral position to mitigate pressure and shear related to total hip procedures.

40110	18 1/4 L x 12 W x 1/2 H"	46.3 x 30.5 x 1.3 cm
--------------	--------------------------	----------------------